


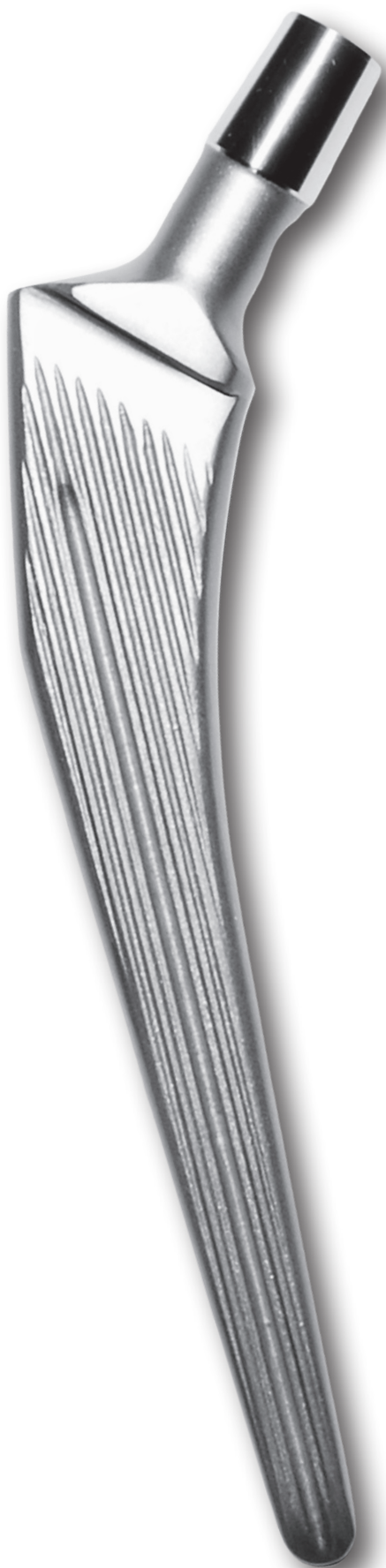

Series 300 / Series 350

Müller Type Straight Stem

Series 300 / Series 350 Options

Stem Design

- Series 300 Müller Type Straight Stem (High-N Stainless Steel)
- Series 350 Müller Type Straight Stem (CoCrMo)
- Standard and lateral Version
- 10 sizes
- Taper 12/14

Müller Acetabular PE Cup

- Various Sizes, available as Standard and SnapFit version
- Material: PE (UHMWPE)

Variable Head Designs

- Modular Metal Heads in Dia. 28/32/36, Necks: S - 3XL
- Bipolar Heads: 19 Sizes from Dia. 42 mm up to Dia. 58 mm

Series 300 / Series 350 Options



Series 20 Bipolar Head and
Müller PE Cup





Metal Head



Series 300 / Series 350 Cemented Stem

Series 300 / Series 350 Cemented Hip Stem

Müller Type Straight Stem

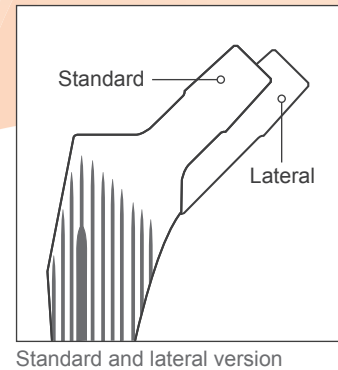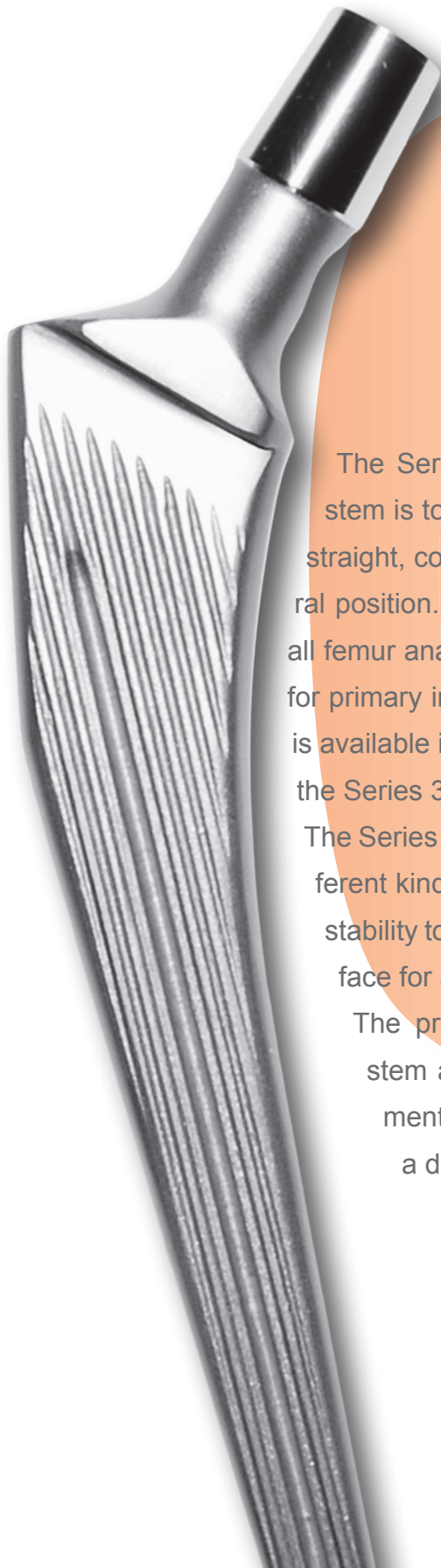

The Series 300 / Series 350 Müller type cemented hip stem is to be inserted without trochanteric osteotomy. The straight, conical shape guides the stem into a neutral femoral position. Due to its flat design, this hip stem suits nearly all femur anatomies either side. The prostheses can be used for primary implantation as well as revisions. The Series 350 is available in 10 sizes, as a standard and lateral version and the Series 300 is available in 11 sizes, as a standard version. The Series 300 / Series 350 Müller Straight stem has two different kinds of grooves. The deep grooves provide reliable stability to the stem. The small grooves create a bigger surface for an optimal mechanical fixation with the cement. The proximal indicated collar prevents a drop of the stem after implantation and compresses the bone cement. The prosthesis consists of a conical design with a doubled elliptical cross section.

Series 300 / Series 350 Stem, Metal Heads & Bipolar Head

Series 300, cemented, sterile

Material: High-N Stainless Steel; Taper: 12/14

Reference Numbers	
Size	Standard
6.25	A1521.05
7.50	A1521.07
8.75	A1521.08
10.00	A1521.10
11.25	A1521.11
12.50	A1521.12
13.75	A1521.13
15.00	A1521.15
16.25	A1521.16
17.50	A1521.17
20.00	A1521.20

Series 350, cemented, sterile

Material: CoCrMo; Taper: 12/14

Reference Numbers		
Size	Standard	Lateral
6.25	A1621.06	A1631.06
7.50	A1621.07	A1631.07
8.75	A1621.08	A1631.08
10.00	A1621.10	A1631.10
11.25	A1621.11	A1631.11
12.50	A1621.12	A1631.12
13.75	A1621.13	A1631.13
15.00	A1621.15	A1631.15
16.25	A1621.16	A1631.16
17.50	A1621.17	A1631.17



Metal Heads, Series 10 Femoral Head

Standard 12/14, sterile

Material: CoCrMo

Neck	Outer Dia.	Ref.
S	28 mm	C1028.01
M	28 mm	C1028.02
L	28 mm	C1028.03
XL	28 mm	C1028.04
2XL	28 mm	C1028.05
3XL	28 mm	C1028.06



S
(-3,5)



M
(+0)



L
(+3,5)



XL
(+7)



2XL
(+10,5)



3XL
(+14)

Neck	Outer Dia.	Ref.
S	32 mm	C1032.01
M	32 mm	C1032.02
L	32 mm	C1032.03
XL	32 mm	C1032.04
2XL	32 mm	C1032.05
3XL	32 mm	C1032.06

Neck	Outer Dia.	Ref.
S	36 mm	C1036.01
M	36 mm	C1036.02
L	36 mm	C1036.03
XL	36 mm	C1036.04

Series 20 Bipolar Head, sterile

Material: Inlay PE, Cup High N

Outer Dia.	Ref.	52 mm	C1252.01
42 mm	C1242.01	53 mm	C1253.01
43 mm	C1243.01	54 mm	C1254.01
44 mm	C1244.01	55 mm	C1255.01
45 mm	C1245.01	56 mm	C1256.01
46 mm	C1246.01	57 mm	C1257.01
47 mm	C1247.01	58 mm	C1258.01
48 mm	C1248.01	59 mm	C1259.01
49 mm	C1249.01	60 mm	C1260.01
50 mm	C1250.01	62 mm	C1262.01



To be used with 28 mm head

Acetabular PE Cup

Müller Acetabular PE Cup

Dia. 28 mm, sterile

Material: PE (UHMWPE)

Outer Dia.	Inner Dia.	Ref.
44 mm	28 mm	B1228.44
46 mm	28 mm	B1228.46
48 mm	28 mm	B1228.48
50 mm	28 mm	B1228.50
52 mm	28 mm	B1228.52
54 mm	28 mm	B1228.54
56 mm	28 mm	B1228.56
58 mm	28 mm	B1228.58
60 mm	28 mm	B1228.60
62 mm	28 mm	B1228.62



Müller Acetabular PE Cup

Dia. 32 mm, sterile

Material: PE (UHMWPE)

Outer Dia.	Inner Dia.	Ref.
44 mm	32 mm	B1232.44
46 mm	32 mm	B1232.46
48 mm	32 mm	B1232.48
50 mm	32 mm	B1232.50
52 mm	32 mm	B1232.52
54 mm	32 mm	B1232.54
56 mm	32 mm	B1232.56
58 mm	32 mm	B1232.58
60 mm	32 mm	B1232.60
62 mm	32 mm	B1232.62

Müller Acetabular PE Cup SnapFit, Dia. 28 mm, sterile

Material: PE (UHMWPE)

Outer Dia.	Inner Dia.	Ref.
44 mm	28 mm	B1328.44
46 mm	28 mm	B1328.46
48 mm	28 mm	B1328.48
50 mm	28 mm	B1328.50
52 mm	28 mm	B1328.52
54 mm	28 mm	B1328.54
56 mm	28 mm	B1328.56
58 mm	28 mm	B1328.58
60 mm	28 mm	B1328.60
62 mm	28 mm	B1328.62



Müller Acetabular PE Cup SnapFit, Dia. 32 mm, sterile

Material: PE (UHMWPE)

Outer Dia.	Inner Dia.	Ref.
44 mm	32 mm	B1332.44
46 mm	32 mm	B1332.46
48 mm	32 mm	B1332.48
50 mm	32 mm	B1332.50
52 mm	32 mm	B1332.52
54 mm	32 mm	B1332.54
56 mm	32 mm	B1332.56
58 mm	32 mm	B1332.58
60 mm	32 mm	B1332.60
62 mm	32 mm	B1332.62

Instruments

Trial Heads, Taper 12/14

Dia. 28 mm	
Neck	Ref.
S	T1728.01
M	T1728.02
L	T1728.03
XL	T1728.04

Dia. 32 mm	
Neck	Ref.
S	T1732.01
M	T1732.02
L	T1732.03
XL	T1732.04

Dia. 36 mm	
Neck	Ref.
S	T1736.01
M	T1736.02
L	T1736.03
XL	T1736.04



Trial Heads for Series 20 Bipolar Head



Outer Dia.	Ref.
42 mm	T1680.42
43 mm	T1680.43
44 mm	T1680.44
45 mm	T1680.45
46 mm	T1680.46
47 mm	T1680.47
48 mm	T1680.48
49 mm	T1680.49


Outer Dia.	Ref.
50 mm	T1680.50
51 mm	T1680.51
52 mm	T1680.52
53 mm	T1680.53
54 mm	T1680.54
55 mm	T1680.55
56 mm	T1680.56
58 mm	T1680.58

Trial Adapter Series 20 Bipolar Head			
S	M	L	XL
T1480.01	T1480.02	T1480.03	T1480.04

Series 300 / Series 350 Mono Rasp

Size	Ref.
7.50	T1241.07
8.75	T1241.08
10.00	T1241.10
11.25	T1241.11
12.50	T1241.12
13.75	T1241.13
15.00	T1241.15
16.25	T1241.16
17.50	T1241.17

Müller Cup Reamers

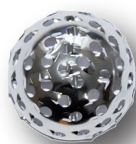

Dia.	Ref.
46 mm	T1340.46
48 mm	T1340.48
50 mm	T1340.50
52 mm	T1340.52

Dia.	Ref.
54 mm	T1340.54
56 mm	T1340.56
58 mm	T1340.58
60 mm	T1340.60
62 mm	T1340.62


Reamer Handle Ref. T1013.01



Stem Impactor/Extractor Ref. T1121.01



Cup Positioner Ref. T1018.01

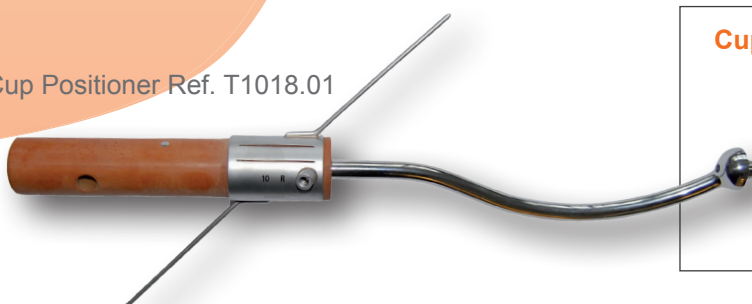

Cup Inserters



Dia.	Ref.
28 mm	T1018.28
32 mm	T1018.32

ARISTOTECH® Implant Precision Forgings GmbH
Im Biotechnologiepark
D-14943 Luckenwalde Germany
phone: +49(0)3371 40640-0
fax: +49(0)3371 40640-215
service@aristotech.de
www.aristotech.de



www.aristotech.de